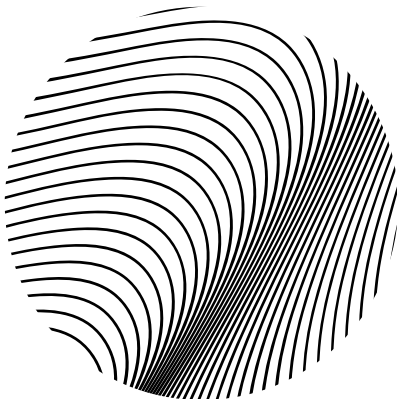
SUPPORTIVE  
TIGHT  
COLLANT  
DE SOUTIEN

# ***HOW TO WEAR THEM***

**Prior to using the device, read the following instructions completely.**

We're stoked you're stoked to get into your new STOKO SUPPORTIVE TIGHTS. Trust us, they're like nothing you've worn before! Here's a couple tips and tricks to make your first experience a little easier.

## **THEY ARE HARD TO GET ON THE FIRST TIME**

Sit down and inch them up like a wetsuit - they will get stuck if you try to pull them on like normal pants. It could take 5 minutes or more to put them on for the first time, so follow the instructions on the next page which walk you through the essential steps to ensure a flawless fit.

## **THEY NEED TO BREAK IN**

They might feel small at first! Stretch the waistband outward and pull on the cables to distribute support down your legs. As you move they will mold perfectly to your body, so the more you move, the more the cables will work their magic.

## **SUPPORTIVE APPAREL FEELS DIFFERENT**

Unlike a traditional brace, your Supportive Apparel will not feel rigid or bulky. The Stoko Support System is designed to mimic your muscles and ligaments, reinforcing your body and engaging when you need it.

## PUTTING ON






1. Pull the dials out to disengage the support system.
2. Stretch the waistband and legs out to relax the cables.
3. Sitting down, pull the tights over each foot. Work the tights up your calf, positioning the crisscrossed cables at your mid-shin. **For all-day boot comfort (such as skiing, hiking, or ice skating), we recommend wearing socks underneath your Supportive Tights.**
4. Continue working the tights up to your knees and align the knee band right below your knee cap.
5. Stand up and pull the waistband over your hips. If the waistband does not feel high enough, pull up on the thigh cables to give yourself some more space. **This will be difficult if the knee band and criss-crossed cables are positioned too low.**





## ADJUSTING

6. Push the dials in and twist to engage support - you will feel them tighten at your waist first.
7. To distribute the support to your legs, pull the waistband outward and move (squat, jog, etc.).
8. STOKO SUPPORTIVE TIGHTS are adjustable and provide more support the more you tighten. Your support is personal, adjust according to comfort or medical recommendation.
9. You will know your tights are working as you feel them tighten across your waist, thighs, and calves.

## REMOVAL

10. Pull the dials out to disengage the support system.
11. Pull the waistband outwards to release the cables and work the tights down your legs and off.

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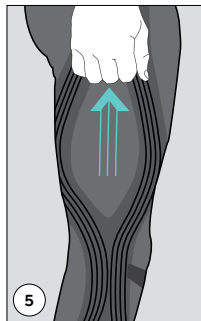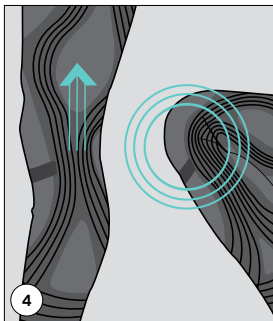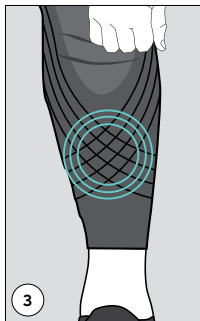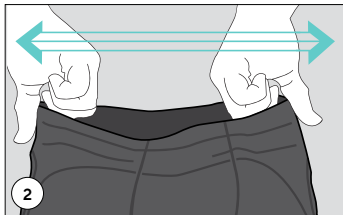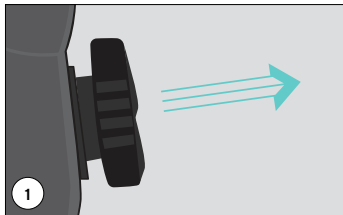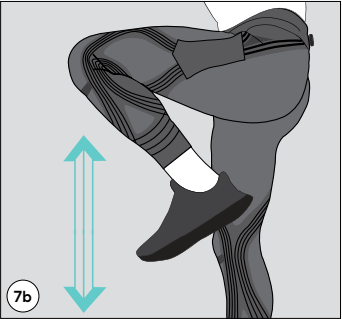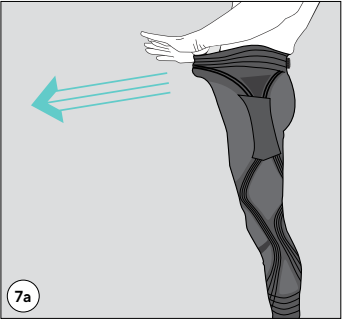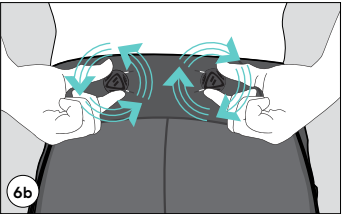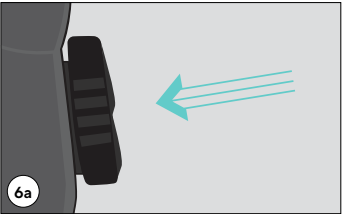# ***THE REST YOU NEED TO KNOW***

## **INTENDED USE**

The STOKO SUPPORTIVE TIGHT is intended for stabilization and support of the knee.

## **INDICATIONS FOR USE**

- Knee pain due to instability
- Knee instability (mild to moderate)
- Knee ligament instabilities (mild to moderate ACL, MCL, and LCL)
- Knee strains and sprains
- Knee osteoarthritis
- Patellofemoral pain syndrome (PFPS)

## **CONTRAINDICATIONS**

If you have any of the following conditions, consult your physician prior to use:

- Pregnancy
- Unstable or displaced fracture of the spine or lower body.
- Conditions exacerbated by compres-

sion or pressure.

- Venous thromboembolism, venous insufficiency or vascular compromise, impaired blood flow, impaired lymphatic drainage, or soft tissue swelling in the lower body
- Impaired sensation of the lower body.

## **WARNINGS**

- Stoko Design Inc. accepts no responsibility for any injury which might occur while using Stoko Design Inc. products or from improper device usage.
- This device is a supportive brace only and is not intended or guaranteed to prevent injury. Consult your health care provider to ensure the STOKO SUPPORTIVE TIGHT is suitable for you.
- Follow all instructions for use of the STOKO SUPPORTIVE TIGHT. Failure to properly size, position, tighten and

adjust the support will compromise performance and comfort. Do not use the STOKO SUPPORTIVE TIGHT over an open wound.

- In rare cases the STOKO SUPPORTIVE TIGHT may cause irritation, discomfort, abrasions, joint pain, or constriction of blood vessels or nerves. In these cases discontinue use immediately. If you experience an unusual reaction consult with a healthcare professional.
- This product does not protect against crushing or hyperextension.
- Do not modify the product. Stoko Design Inc. accepts no responsibility for the consequences of any adjustments, modifications and/or repairs made to the STOKO SUPPORTIVE TIGHT by any user or other third party.

## **DUTY TO REPORT**

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

## **INTENDED FOR SINGLE PATIENT USE.**

### **USE LIFE**

Discontinue device use after 5 years or if you see significant signs of wear or failure.

### **PRODUCT CARE**

Keep your STOKO SUPPORTIVE TIGHT fresh! Follow the instructions on the sewn-in-label, wash in the washbag provided, and hang to dry. Do not expose the product to direct heat.

### **MATERIAL CONTENT**

Garment: 64% Nylon / 36% Elastane  
Support System: UHMWPE, Nylon, Polycarbonate, Polyoxymethylene, Polyurethane.

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